

PRODUCT BROCHURE

KERN SET-18

Date: 19.04.2025

ENGLISH

Software EasyTouch

SET-18 Container Master

EasyTouch Container Master – Centralised container management

CONTACT






KERN & SOHN GmbH
Ziegelei 1
72336 Balingen
Germany

Phone : +49 7433 9933-0
Fax : +49 7433 9933-149
Email : info@kern-sohn.com
Web : www.kern-sohn.com



Scan here for more
information

PRODUCT SHOWCASE

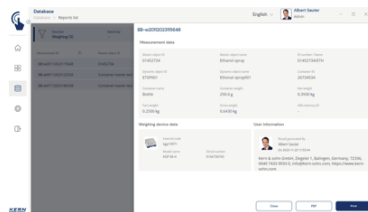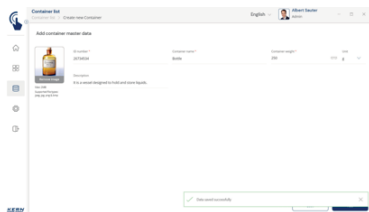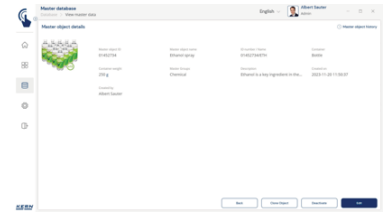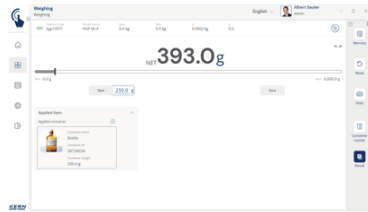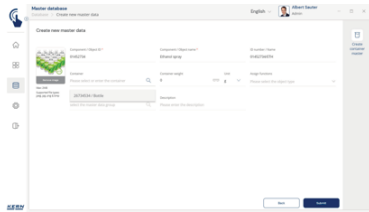

DESCRIPTION

- Note: Prerequisite for this set is the basic program SET-01 Base
- Container management: Different items are stored and weighed in similar containers. The tare weights of these containers can be managed centrally with SET-18 Container Master. If the weight of a container changes, the new tare weight is simply updated in Container Master and automatically applied to all associated items and functions in the application
- Fixed assignment: A container can be permanently assigned to each master data object to ensure consistency
- Container display: When using a master data object, in EasyTouch users see the permanently assigned container or a manually selected container. This reduces the risk of confusion and improves accuracy in the work process
- Container selection: Users can also select a container during the weighing process without having to define it beforehand, for example by scanning the container ID with the barcode. This helps to avoid incorrect tare entries

SPECIFICATIONS





Article codes, product groups & Status information

Model article-code	SET-18
Model-Serie	ETS
Assortment	KERN
Product Type	Software
Product group	EasyTouch
Product usage type	Main item;Accessories/Services
EAN code (Model)	404576 1478312
Customs tariff number	852349 10

Technical data - Packing & Shipping

Delivery time	1 d
---------------	-----

FUNCTIONS

Standard
   
Option
